



TECHNIMOUNT SYSTEM™



BRACKET PRO SERIE® 80-FL

MOUNTING SYSTEM FOR AIR AND GROUND TRANSPORT OF THE
B.BRAUN INFUSOMAT® VOLUMETRIC AND PERFUSOR® SPACE PUMPS

SAFETY & FLEXIBILITY WHERE IT MATTERS MOST

The Bracket Pro Serie® 80-FL is designed for Air and Ground EMS and Critical Care Transport of the B.Braun Infusomat Volumetric and Perfusor Space Pumps.

Lightweight and robust, it can secure stacks of 3 or 4 B.Braun pumps in different applications including surface, stretcher/cot and wall systems.



Ambulance



Stretcher



Military



Air

At Technimount, we understand that EMS Providers need solutions that they can rely on. We design and develop quality mounting systems which are safe and flexible, allowing you to concentrate on what matters most - patient care.


BRACKET PRO SERIE® 80-FL

POTENTIAL CONFIGURATIONS

Safety MD Transporter



Wall Mount System (EU)

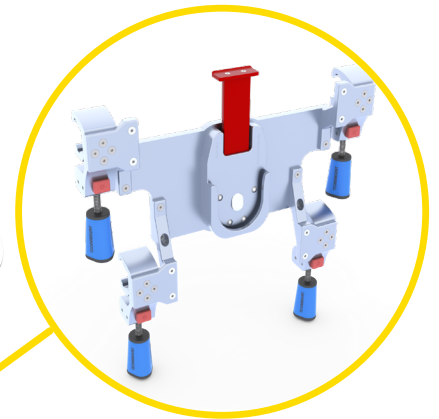

Safety MD Transporter
(Horizontal)



Proplate with 2 Micro
Bases for IT Stretchers

TECHNICAL SPECIFICATIONS

Product Name	Bracket Pro Serie® 80-FL
Compatibility of Medical Device	B.Braun Infusomat® Volumetric and Perfusor® Space Pumps
Environment of Usage	EMS/CCT (Ground & Air) – Hospital/Clinic, Military
Device Section Material	<ul style="list-style-type: none">• Bracket - High-Density Aluminum with Blue-Anodized Finish• Shelves - UV-Stabilized White ABS• IV Pole - High-Density Aluminum with Silver-Anodized Finish• IV Pole Hooks - Stainless Steel
Bottom Micro Disc Section Material	High-Density Aluminum with Silver-Anodized Finish
Anti-Rotation Micro Disc Material	High-Density Aluminum with Silver-Anodized Finish
Temperature of Usage	- 35° C @ 45° C (Direct Conduction)
Impact Resistance – Ground	SAE J3043 (30 G) Plus 10% Safety Factor (in process)
Impact Resistance - Flight	FAA – FAR -DOT 27.561 & DOT 27.562 up to 20G (in process)
Installation	Proprietary Locking Mechanism on Each Shelf
Attachment Options	2 Locking Discs – Bottom Micro Disc – Anti-Rotation Micro Disc
Transportation Compatibility System	Compatible with Technimount's: <ul style="list-style-type: none">• Wall Mount Systems• Safety MD Transporters• Proplate with Micro Bases for Incubators
Options	<ul style="list-style-type: none">• Bracket Pro Serie 80-FL for Stacks of 3 B.Braun Pumps• Bracket Pro Serie 80-FL for Stacks of 4 B.Braun Pumps
Features	<ul style="list-style-type: none">• Robust Mounting Solution for EMS/CCT/ALS• Air and Ground Compatibility• Multiple Configuration Options for Transport• Red Quick Release Mechanism to Easily Insert and Remove the Pumps• Polyurethane Plastic Transport Handle
Physical Overall Dimensions	Bracket Pro Serie 80-FL3 <ul style="list-style-type: none">• Height 13.40 in., Width 10.45 in., Depth 8.40 in. Bracket Pro Serie 80-FL4 <ul style="list-style-type: none">• Height 16.48 in., Width 10.45 in., Depth 8.40 in.
Weight	Bracket Pro Serie 80-FL3 <ul style="list-style-type: none">• Bracket - 8.58 lbs• Bracket with 3 B.Braun Pumps - 17.82 lbs Bracket Pro Serie 80-FL4 <ul style="list-style-type: none">• Bracket - 11.52• Bracket with 4 B.Braun Pumps - 23.84
Part Number	<ul style="list-style-type: none">• Bracket Pro Serie 80-FL3 - 1110-11-BBRK3FL• Bracket Pro Serie 80-FL4 - 1110-11-BBRK4FL• Micro Base - 120-20-UN

YOUR SAFETY IS OUR PRIORITY



FEATURES

- ✓ Designed to secure 3 or 4 B. Braun Infusomat or Perfusor Space Pumps during Air and Ground EMS/Critical Care Transport
- ✓ Medical device retained on main assembly with proprietary locking mechanism on each shelf
- ✓ Red Quick Release Mechanism to instantly remove the pump from the bracket
- ✓ Anti-Rotation Micro Disc and Bottom Micro Disc to support multiple mounting configurations on surface, stretcher/cots and wall systems
- ✓ IV Pole with Quick Release Mechanism for use during temporary transfer from bed to stretcher
- ✓ High-density aluminum with blue-anodized finish
- ✓ Designed to comply with FAA - FAR - DOT 27.561 & DOT 27.562 up to 20G (in process) and SAE J3043 (certification in process)

BENEFITS

- + Flexible for use in different applications (surfaces, stretcher/cot and wall systems) and adaptable to quantity of pumps required for transport
- + Easily install and remove the medical device from the bracket for maintenance
- + No interference with the medical device components, cables and tubing
- + Access to medical device from all angles (360° rotation on the Safety MD Transporter with Vertical and Horizontal Micro Base)
- + Quickly install the IV Pole during temporary transfer from bed to stretcher
- + Certified mounting system that complies with the highest industry standards and regulations
- + Reduces risks of damage to the equipment and ensures the safety of patients and crew during EMS and Critical Care transport



Bracket Pro Serie 80-FL4 for Stacks of 4 B-Braun Pumps

SAFETY AND FLEXIBILITY WHERE IT MATTERS MOST

technimount.com

